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BRESSER Science ADL-601F
Transmission-type and
Epifluorescent Microscope

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CAUTION!

To work with this microscope, sharp and pointed aids are being used. Please take care that this microscope and its accessories are stored at a place out of reach of children. Let children only work with this microscope under an adult's supervision! Keep packing material (plastic bags etc.) away from children!

Note!

If you have any complaints or queries please first contact your national service centre by telephone. The address is included in these instructions.

Only for EU countries

Do not dispose of electric equipment together with household waste material!

In observance of European Directive 2002/96/EC on waste electrical and electronic equipment (WEEE) and its implementation in accordance with national law, electric equipment that have reached the end of their life must be collected separately and returned to an environmentally compatible recycling facility.

Discharged batteries and damaged re-chargable batteries must be disposed of at special battery collection points. Information is available from your local disposal agent or local authority regarding the disposal of devices or batteries manufactured after the 01.06.2006

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WARNING


Please disconnect the plug from mains socket before changing the fuse, removing the bottom or the rear plate or opening the lamp door to change the lamp. The power cord provided with the equipment has a grounded plug. Always use the power cord with a properly grounded mains socket.

Do not expose the instrument to high temperatures or humidity. Avoid using the instrument in extremely dusty locations.

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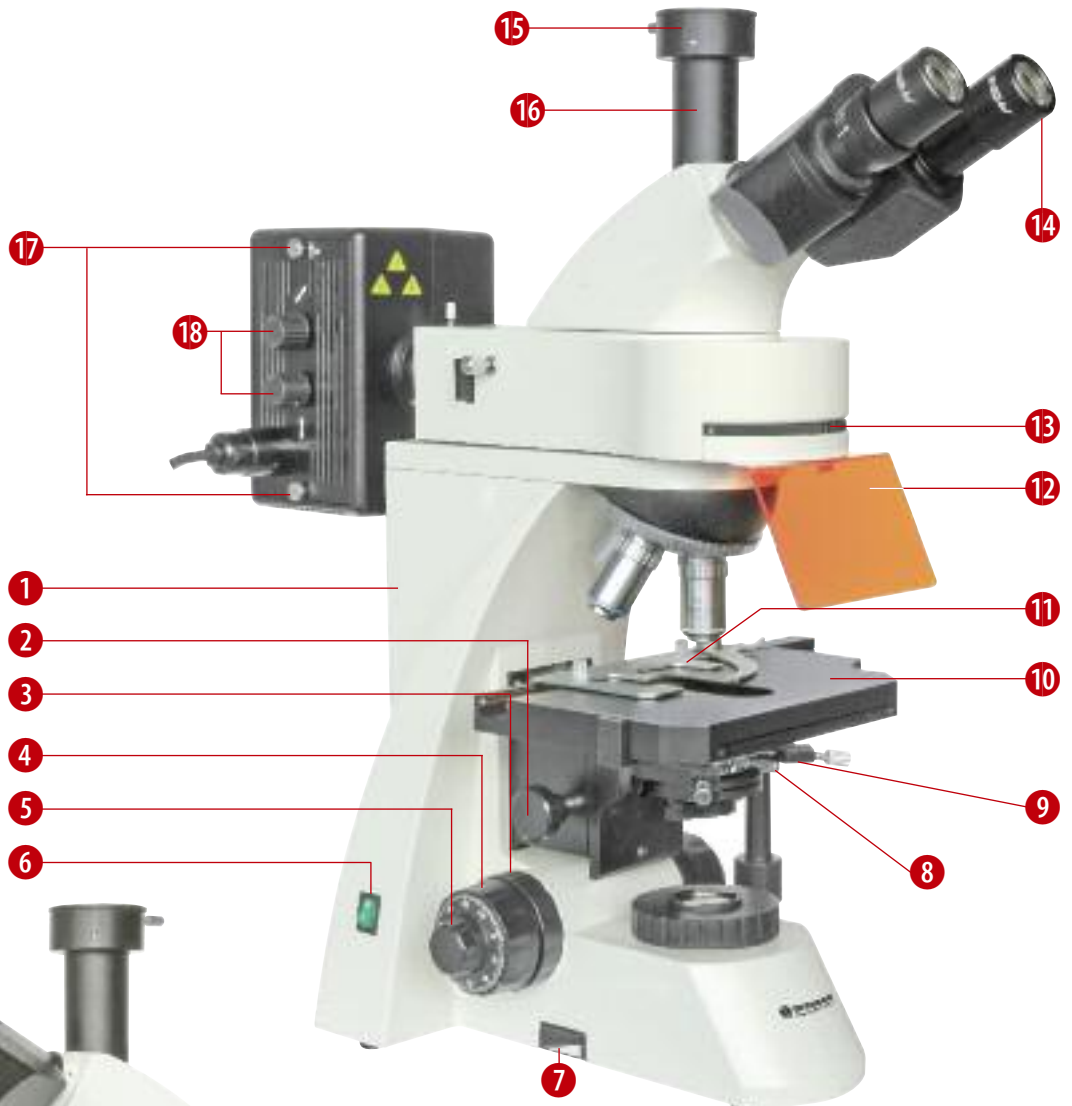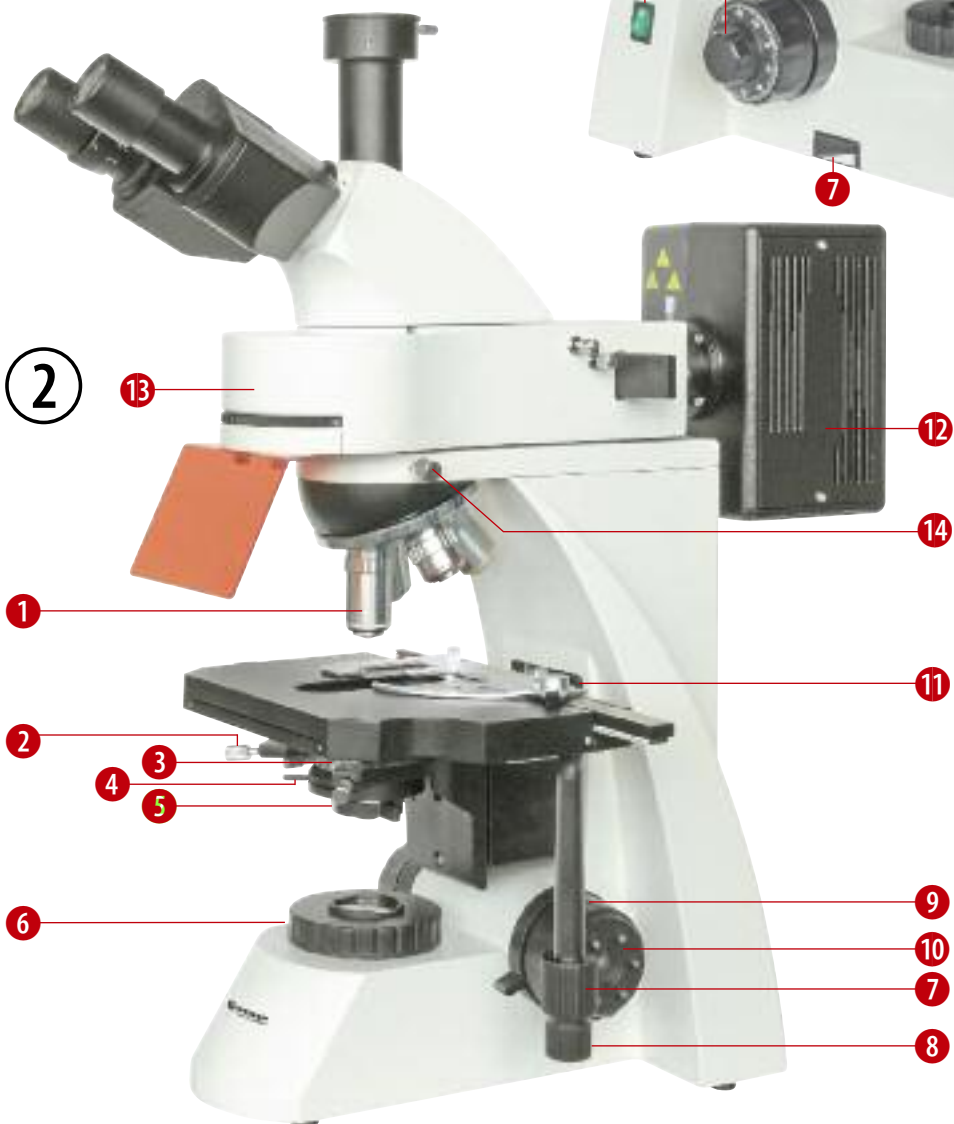
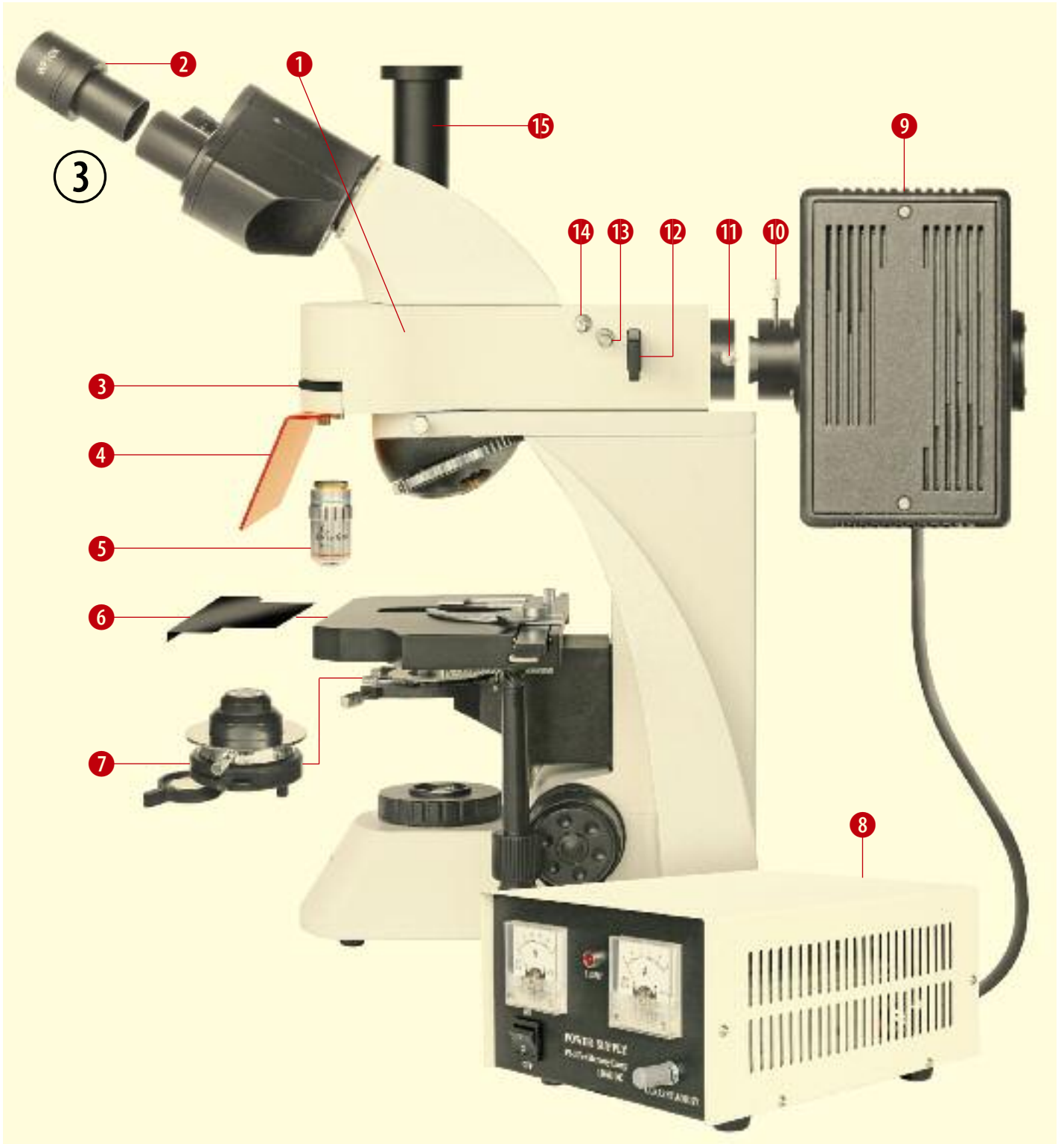
Table of Contents

- I. Components
- II. Technical Data
- III. Installation
- IV. Observing Operation
- V. The Operations of Each Unit
- VI. Maintenance, EEC Conformity Explanation, Warranty

I. COMPONENTS

Fig. 1

- 1 Stand
- 2 Condenser up-down knob
- 3 Knob for tensional adjustment of focusing
- 4 Coarse focusing knob
- 5 Fine focusing knob
- 6 Transmitted light power switch
- 7 Transmitted light dimmer
- 8 Condenser holding screw
- 9 Condenser holder
- 10 Stage
- 11 Slide holder
- 12 Radiation-proof baffle
- 13 Filter wheel
- 14 Eyepiece
- 15 Camera port
- 16 Photo tube
- 17 Side-board retaining screw of the lamp box
- 18 Lamp centering knob

Fig. 2

- 1 Objective
- 2 Condenser centering screw
- 3 Condenser mounting ring
- 4 Aperture diaphragm adjustment lever
- 5 Auxiliary condenser lens
- 6 Field diaphragm
- 7 Cross moving knob
- 8 Lengthwise moving knob
- 9 Fine focusing knob
- 10 Coarse focusing knob
- 11 Limit knob (up stop)
- 12 Mercury lamp box
- 13 Microscope head with incident-type fluorescence illumination unit
- 14 Microscope head holding screw

II. TECHNICAL DATA

The "BRESSER Science ADL-601F" is used for incident-type fluorescence microscopy and for transmission-type bright field microscopy. By using incident light excitation and observing binocularly, the microscope provides a clear and high contrast fluorescent image. It is an ideal instrument for

Fig. 3

- 1 Microscope head (with incident illumination unit): Take off the protective cap under the head. Loosen the microscope head holding screw (on the right side above the nose-piece) and install the head to the frame by inserting the dovetailed flange (on the underside of the head) into the frame. Secure the head by tightening the holding screw. Do not release the head from your grasp until it is firmly secured to the frame!
- 2 Eyepiece: Take off the protective caps from the eyepiece holders and insert both eyepieces into the eyepiece holders.
- 3 Filter wheel: For transmission-type bright field microscopy, set it to the position "0"; for incident-type fluorescence microscopy, set it to the position of the chosen filter (UV, V, B, G).
- 4 Radiation-proof baffle: Screw the baffle on the front of the head.
- 5 Objective: Install the objectives in the nosepiece.
- 6 Light stop plate: This plate can be pushed into the front of the stage.
- 7 Condenser (with auxiliary condenser lens): Install the condenser in the condenser holder.
- 8 Mercury lamp power supply unit: Plug in the power supply plug for the mercury lamp into the power supply socket on the rear of the unit and tighten the knurled ring. Plug in the power supply plug of the power cord into the power supply socket of the mercury lamp power supply unit on the rear of the unit.
- 9 Mercury lamp box: With the holding screw (10), install the lamp box on the incident illumination unit. Install the mercury lamp in the lamp box.
- 10 Condenser adjustment lever.
- 11 Holding screw of mercury lamp box.
- 12 Filter rod.
- 13 Field diaphragm centering screw.
- 14 Field diaphragm adjustment rod.
- 15 Photo tube. For photo attachment (optional): Take off the protective cap from the photo tube and install the photo attachment.

After checking and making sure that the voltage of the instrument is in accord with the supplied voltage, plug the mains plug into the mains socket.

research of biology, cytology, oncology, genetics, immunology etc. At the same time, it is also applied for ranges such as analysis of sedimentary rock, inspection of impurity of semiconductor, environmental protection, microchemistry etc.

NOTE

The mercury lamp in the mercury lamp box had been taken down and placed in a foamy packing box in order to prevent it from damage in the period of transportation. Please replace the stick in the mercury lamp box by the mercury lamp. (Refer to chapter "V. OPERATIONS OF EACH UNIT – Exchanging the mercury lamp for fluorescence illumination".) Please store the stick for further use.

TECHNICAL DATA

Eyepiece

Eyepiece type	Magnification	Field of view Ø / mm	Remarks
Wide field eyepiece (WF)	10x	20	
WF eyepiece with crosshair	10x		optional

Objektive

Objective type	Magnification	Numerical aperture NA	Remarks
Fluorescence objective	40x (glycerin)	1	
Planachromatic objectives	4x	0.1	
	10x	0.25	
	40x	0.65	
	100x (oil)	1.25	for oil immersion (preparation without cover glass)

Total magnification

Objectives:	4x	10x	40x	100x
Total magnification				
Eyepiece:				
10x	40x	100x	400x	1000x

- Numerical aperture of condenser: 1.25.
- Moving range of mechanical stage: 50 mm x 75 mm.
- Minimum division of fine focusing scale: 0.002 mm.
- Adjustment range of interpupillary distance: 53 mm – 75 mm.
- Incident light source for fluorescence: 100 W super-high pressure mercury lamp with DC power supply (AC 230 V / 50 Hz).
- Transmitted light source: 6 V / 20 W halogen bulb lamp with dimmer and power supply (AC 230 V / 50 Hz).
- Exciting filters: Ultraviolet (UV), Violet (V), Blue (B), Green (G).

III. INSTALLATION


Remove all parts from their packing materials and retain the packaging in the event you need to transport the product. Install all the units according to Fig. 3.

IV. OBSERVING OPERATION

1. Transmission-type bright field microscopy (Fig. 1)

1. Press the power switch (6) of "I", which means that power supply for transmitted light illumination is switched on.
2. Put a preparation on the stage (10), turn the filters turnplate (13) to position "0" (clear opening) and the 10x objective (Fig. 2, No. 1) to the working position and bring the preparation into exact focus using the coarse and fine focusing knobs (Fig. 2, No. 9&10).
3. Adjust the interpupillary distance and the diopter when you are watching binocularly with both eyepieces (14).
4. Adjust the height position of the condenser (Fig. 3, No. 7) using the condenser up-down knob (2). Adjust the transmitted light dimmer (7) and the aperture diaphragm on the condenser (Fig. 2, No. 4) in order to get a satisfying illumination. If you watch with the 4x objective, turn the auxiliary condenser lens (Fig. 2, No. 5) into the optical path so that you can get symmetrical light.
5. Adjust the fine focusing knob (5) when you change between various objectives to observe.
6. Place some immersion oil between the 100x objective and the cover glass and between the condenser and the slide, when you observe with the 100x objective.

2. Incident-type fluorescence microscopy (Fig. 4)

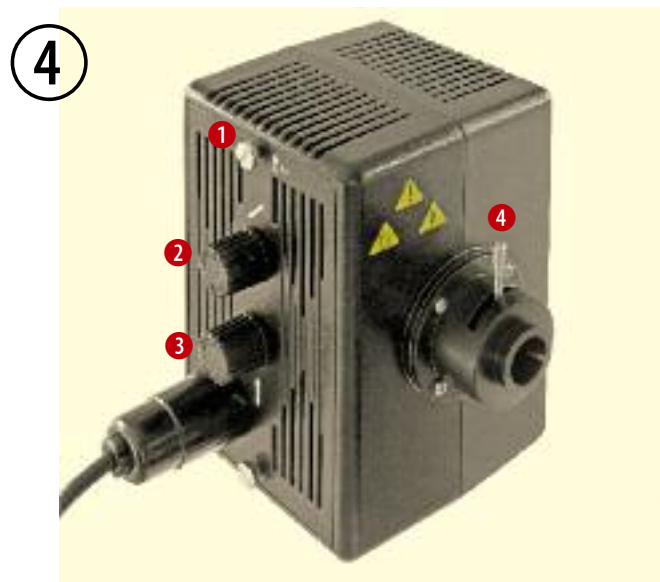1. Turn the filters turnplate (Fig. 1, No. 13) to the position corresponding to the chosen excitation radiation (UV, V, B or G, not 0).
2. Switch on the power supply unit (Fig. 3, No. 8) by pressing the power switch to "I/ON". If the indicator "LIGHT" is bright, it shows that the mercury lamp is lighted up.
3. Place a piece of white paper on the stage and remove the objective on the working position. Adjust the condenser adjusting lever (4) to image the mercury lamp arc on the paper without objective.

4. If the image of the arc is not centered in the overpass hole, turn the lamp left/right centering knob (2) and the lamp up/down centering knob (3) on the lamp box to ensure that the image of the arc is in the center of the field.
5. Install the removed objective again, interchange the 10x objective into the working position and use it to observe. Adjust the coarse and fine focusing knobs (Fig. 1, No. 4&5) to image the field diaphragm clearly. Use the field diaphragm adjustment rod (Fig. 3, No. 14) to close the field diaphragm small, and then adjust the field diaphragm centering screws (Fig. 3, No. 13, on both sides) to center the field diaphragm. Use the field diaphragm adjustment rod again to make the image of the field diaphragm slightly larger than the field of view in the eyepiece. The 40x fluorescence objective may be used for fluorescence microscopy.
6. If the field appears too bright, pull the matted glass filter on the filter rod (Fig. 3, No. 12) into the working position.


NOTES:

- a) For use of the mercury lamp, always let be an exciting filter (UV, V, B or G) in the optical path, make sure that the whole housing of the instrument is closed, and do not look into the eyepieces when you change between the excitation filters.
- b) After lighting up, the mercury lamp needs a minimum of 15 minutes to work stably.
- c) Wait at least 10 minutes for bulb cooling after turning off, then the lamp can be restarted.
- d) Work with glycerin immersion when you use the 40x fluorescence objective.
- e) Pull the filter rod (Fig. 3, No. 12) blocking into the optical path to protect the object from light and radiation and to avoid damaging an observation result.
- f) After turning off the mercury lamp, turning the filters turn-plate (Fig. 1, No. 13) to the clear opening position "0" and removing the light stop plate (Fig. 3, No. 6), transmission-type bright field microscopy can be done.
- g) Notes to some risks at the use of mercury high pressure lamps:

You can find notes to some risks at the use of mercury high pressure lamps on page 19.

V. THE OPERATIONS OF EACH UNIT

1. Interpupillary distance adjustment (Fig. 5)

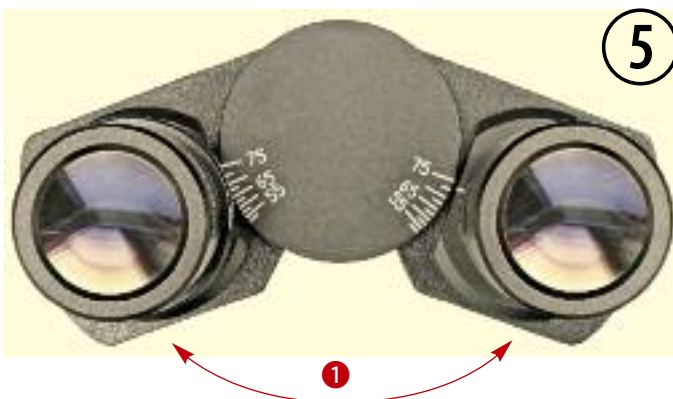Put a specimen on the stage and focus to obtain a clear image of it. Adjust the interpupillary distance of the binocular tube by "folding" (1) until the right and left field of view can be composed to one.

2. Diopter adjustment (Fig. 6)

Put a specimen on the stage. Turn the 40x objective to the working position. Firstly, observe through the right eyepiece with the right eye; adjust the coarse and fine focusing knobs to image clearly. Secondly, observe through the left eyepiece with the left eye; adjust the diopter adjustment ring (1) to image clearly.

3. Coarse/fine focusing (Fig. 7)


The instrument uses a coaxial coarse/fine focusing mechanism to raise or lower the stage. (2) is the coarse focusing knob for coarse adjustment, (1) is the fine focusing knob for fine adjustment. The smallest graduation on the scale of the fine focusing knob is 2 μm of vertical.

The tension of the coarse focusing knob is adjustable and pre-set at the factory for ease of use and to prevent the stage from naturally sliding down. If you wish to adjust the tension, use the knob for tensional adjustment of focusing (3). Turning the ring toward the front of the microscope increases the tension, and toward the rear of the microscope loosens it. Tension is too high if you experience physical discomfort.


The limit knob (4) prevents accidental contact between longer objectives and the specimen or the stage. After locking the limit knob (by turning up and fastening) at a chosen height of the stage, you cannot focus the object closer to the objective using the coarse focusing knob. So, the object is protected against damage. Its use also simplifies focusing (so-called pre-focusing). After changing specimens or objectives, focusing is easily accomplished by rotating the coarse focusing knob to reach the pre-focused position. Then make fine adjustments with the fine focusing knob. Focusing movement with the fine focusing knob isn't affected by using the limit knob.

4. Stage (Fig. 8)

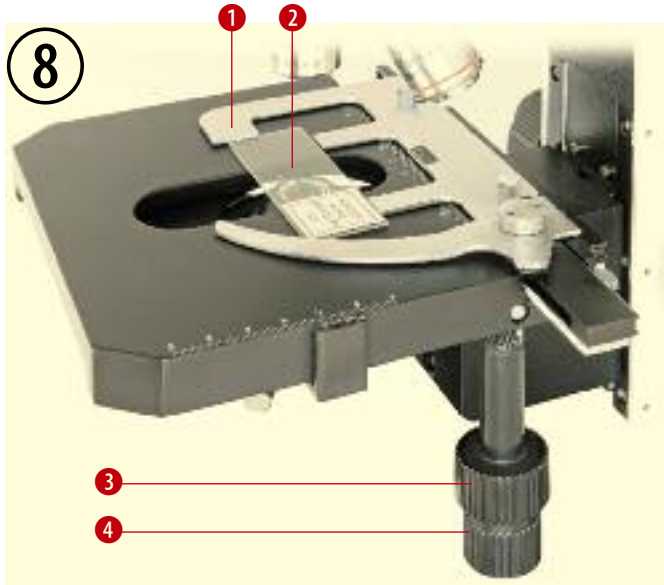The convenient slide holder (1) on the stage is used for holding the slide glass (2). The longitudinal (lengthwise) adjustment knob (3) and the transversal (cross) adjustment knob (4) are coaxial. The stage moves expediently.

5. Adjustable condenser (Fig. 9)

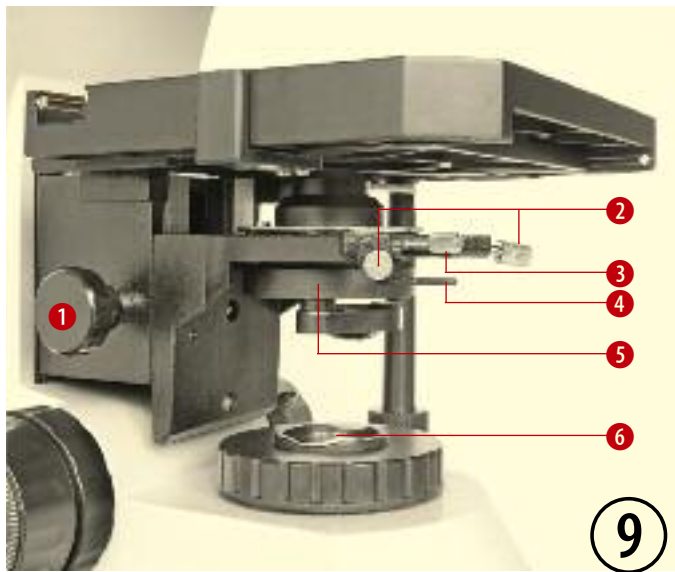The condenser (5) is moved up or down via turning the up-down knob (1). For centering, the condenser centering screws (2) are used. The condenser can be taken down easily if one unscrews the condenser holding screw (3). The place for a filter plate is on the filter seat on the field diaphragm (6).


6. Power switch and adjustable brightness (Fig. 1)

Turn on the power switch (6). Adjust the dimmer (7) until the image can be observed comfortably. Note: Don't let the dimmer at the position of maximum brightness for a long time. That reduces the life-span of the lamp.

7. Aperture diaphragm (Fig. 9)

The aperture diaphragm adjustment lever (4) can be turned in order to open or close the aperture diaphragm. Remove an eyepiece and watch through the eyepiece tube. You will see the bright aperture diaphragm image in the dark objective pupil. For centering, the condenser centering screws (2) are used when the aperture diaphragm image is eccentric with the objective pupil. Turn the aperture diaphragm adjustment lever (4) for getting a good resolution and contrast perception. Usually, the diameter of the aperture diaphragm image, which has to be adjusted, is 70-80 percent of the objective pupil.

8. Exchanging the mercury lamp for fluorescence illumination (Fig. 10)

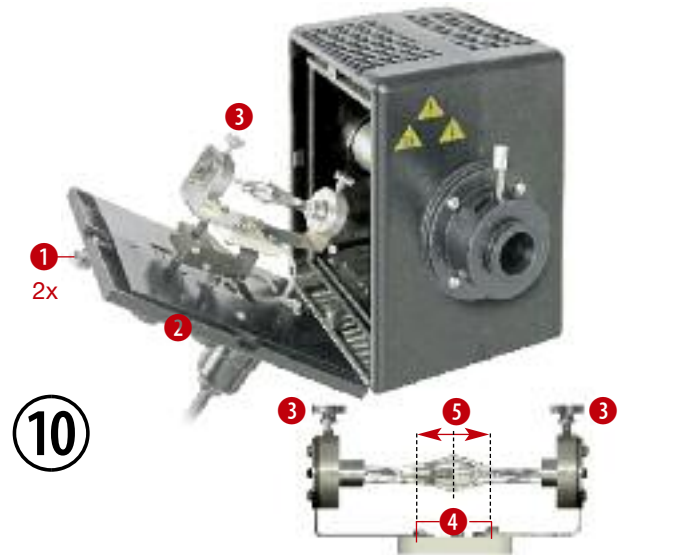1. Switch off the mercury lamp power supply unit (Fig. 3, No. 8), unplug the mains plug. The mercury lamp has to be cooled down.
2. Loosen the sideboard retaining screws of the lamp box (1) and remove the sideboard (2).
3. Loosen the two lamp holding screws (3) located at the both ends of the lamp. Put down the mercury lamp (5).
4. Clean the new lamp with anhydrous alcohol. Avoid fingerprints on the lamp, don't touch the lamp with bare fingers!
5. Assemble the new mercury lamp (5). Its center should be located over the two little screws (4). Tighten the lamp holding screws (3).
6. Mount the sideboard (2). The center of the lamp can be moved with the two lamp centering knobs (Fig. 4, No. 2&3).


9. Exchanging the fuse for fluorescence illumination (Fig.11)

1. Switch off the mercury lamp power supply unit (Fig. 3, No. 8) and unplug the mains plug of the unit to disconnect the power supply.

- On the rear of the power supply unit and next to the mains in, there is the screw of the fuse holder. Loosen the screw of the fuse holder (push and turn left), fetch out the bad fuse, replace a new fuse (230 V / 5 A, F) and tighten up the screw of the fuse holder.


10. Exchanging the halogen lamp for transmitted light illumination (Fig. 12)

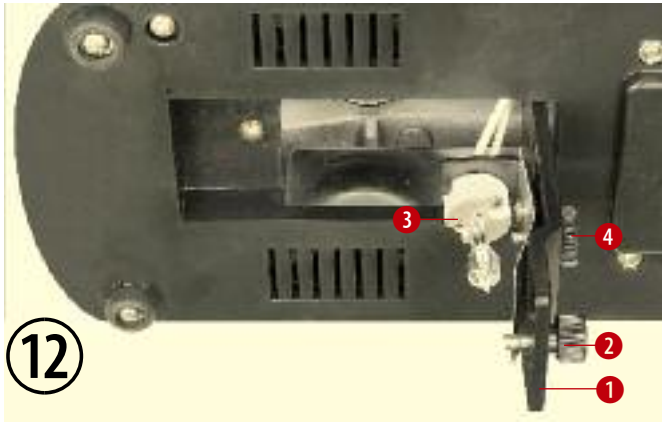- Switch off the transmitted light power switch (Fig. 1, No. 6) and unplug the mains plug of the microscope to disconnect the power supply. The halogen lamp has to be cooled down.
- Incline the microscope, loose the lamp door retaining screw (2) on the lamp door (1) and open the lamp door.
- Pull out the old lamp from the lamp base (3).
- Clean the new lamp with anhydrous alcohol; don't touch the lamp with bare fingers.
- Put the new lamp into the lamp base (3). Notice the properly touching.
- Reinstall the lamp door (1) on chassis with the lamp door retaining screw (2).
- Mount the lamp well, plug in the mains plug, turn an objective into the optical path, adjust the condenser up and down, and make light enter the view field. If the illumination is non-uniform, loose the lamp centering screw (4) slightly and move the lamp base (3) with it to bring the lamp spot into the center. Finally, tighten up the screw (4).

11. Exchanging the fuse for transmitted light illumination (Fig. 13)

- Switch off the transmitted light power switch (Fig. 1, No. 6) and pull out the plugs of the power cord from mains socket and from mains in at the microscope (1).
- Grip with a small screwdriver under the lug of the fuse holder (2) and pull it carefully forward. Take out the fuse holder (2). Change the fuse. The specification of the fuse: 250 V / 0.5 A, F. Install the fuse holder and plug again.

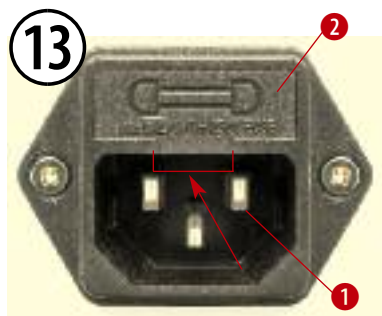VI. MAINTENANCE

1. Cleaning the frame and the stage

Disconnect the plug from mains socket before cleaning. Clean the frame and the stage with a soft cloth moistened with a mild detergent solution. Be sure that the instrument is dry before using.

2. Cleaning optical parts

The eyepieces and the objectives of the microscope are coated. They should not be wiped because dry dirt or dust may scratch the coating. It is best to remove the parts to be cleaned from the frame prior to cleaning. Always blow loose dust away first. Use lens tissue of good quality moistened with a lens cleaner or a small amount of alcohol; then wipe the surface clean with a lens tissue. Solvents such as xylene should NOT be used as cleaner.

3. Do not disassemble objective lenses!

4. Cleaning the 100x oil immersion lens

The immersion oil should be removed from the lens at the end of each workday using a lens tissue moistened with a lens cleaner or a small amount of alcohol.

5. Being not used

Cover the microscope with the dust protective cover and place it there, where it is dry and mildewless. We suggest the storage of all objectives and eyepieces in a closed container with drying agent.

Risk	Possible effect	Counteraction
Power supply of the lamp	Short circuit, lamp box energized	a) Never open or touch an energized lamp box b) Plug off the mains plug, separate all power supply from mains (also do so for installations at the lamp box or at the lamp)
Lamp is very fragile	Break of the lamp	Install the lamp very careful without pressure
Lamp produces a lot of heat and becomes very hot	Danger of burn (your fingers etc.)	Before opening, let the lamp box cool down minimum 15 minutes
	Danger of combustion	Always keep open the aeration slots; don't cover them
Lamp produces light and ultraviolet radiation (UV) of high intensity	Glare; permanent (but not really coming up at once) damage of the eyes (cornea, lens, retina) up to blindness; burn of the skin	a) Shielding from light and radiation (e. g. by a closed lamp box and, in general, a closed housing; use of a blocking filter, a light stop plate and an orange radiation-proof baffle); never open the lamp box when the lamp is lighting! b) Personal safety (never look into the lighting lamp, the leaving ray (neither with protective goggles), into its reflection or image; if necessary, wear UV protective goggles and gloves)
UV radiation of the lamp produces poisonous gases in air (ozone, oxides of nitrogen)	Already at relative low concentration of gas: "smell of ozone", dryness in the nose, feeling of burn in the throat, headache, nausea, irritations of the mucous membranes	Operate the lamp in a large room with good ventilation
Overheating of the lamp	Danger of explosion	a) Make sure the correct power supply of the lamp b) Adjust the lamp immediately after installation
Lamp reaches the end of the lifetime (discoloring of the lamp)	Increasing danger of explosion	Exchange the lamp early enough
Explosion of the lamp	Flying splinters, emission of hot gas	Shielding of the lamp by the closed lamp box
	Emission of mercury, which can be inhaled or received by skin and cause damage to someone's health	a) Leave the room immediately to avoid inhalation b) Air the room carefully (minimum 20 – 30 min) c) After cooling down: Remove thoroughly the deposited rests of mercury (e. g. by commercially available adsorbent agents)
Lamp contains mercury	Environmental damages by the mercury	Disposal of the old lamp only in accordance with the national or European prescriptions (in Europe: as hazardous waste under EWC code 20 01 21* "Fluorescent tubes and other mercury containing waste")
Transport of the instrument	Break of the lamp, damage in the lamp box	a) Remove the lamp and transport it in the original packing box b) Reinstall the stick (transportation lock) in the lamp box to protect the moveable parts there

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Bresser GmbH